

Mobile, adjustable and self contained workstations designed to enhance lab productivity. Push a button and the work surface adjusts to meet the height demands of individual users and equipment. Add plumbing, shelving and accessories to meet your specific requirements. You can easily reconfigure lab layouts as needs change; simply move it where you need it when you need to.



OPTIMA™

MOBILE WORKSTATION

Highly Adaptable - Can be conveniently integrated into existing work spaces or located as a self supporting workstation.

Integrated Services - Plumbing, electrical and data services are factory installed and designed to quick connect to any of Mott's overhead service carriers.

Upper Storage - Optional fixed shelves mount between the uprights to optimize vertical space and allow for maximum utilization of valuable work surface space.


Removable End Panels - Designed for rugged durability, these panels provide an inspection point for connections, shut-offs and service lines. Access for modifications is quick and easy.

Service Accessibility - Hinged panel is easily accessed for convenient service maintenance.

Flexible Solutions - Easily connects to overhead service carriers as needed and moved out of the way when finished.

Configurations - Available in single sided and double sided assemblies in various widths.

Accessories - Optima™ Workstations can be configured with many accessories such as keyboard trays, monitor arms, and shelf task lights to suit your needs.

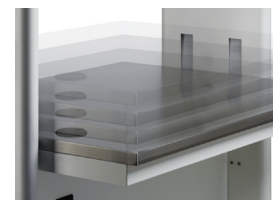

Adjustable Work Surface - Electrohydraulic height adjustable table frame with motor has the flexibility to be raised or lowered from 30" to 36" at the push of a button to accommodate the requirements of individual users, different procedures or equipment.



Plumbing Services - Pre-plumbed service fixtures are a standard option to supply the necessary utilities and are easily accessible while keeping the work surface open.



Effortless Mobility - Swivel casters have an adjustment knob that lowers the pad and raises the wheel off the ground to prevent rolling. Perfect for applications where work demands stability and a level work surface.

mott
MANUFACTURING

ISO 9001:2008 REGISTERED

QUALITY BY DESIGN